

**2011
NEW!!**

ATTO Blocking solution series ! EzBlock Chemi / BSA / CAS

- Higher efficiency of blocking than conventional product !
- Various line-up !
- Easy preparation in 5-fold dilution with distilled water !
- Adjustable of Tween-20 according to purpose of use !


Preparation

- Dilute in 5-fold with distilled water

For EzBlock BSA / CAS, mix an appropriate volume of Tween-20 included as a standard composition.

(ex) Use 50mL blocking solution for 10cm x 10cm membrane.

Blocking Data of New Blocking Solution

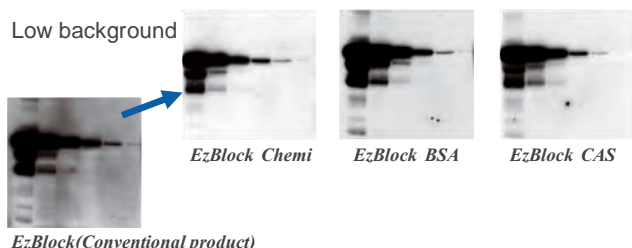

Features

EzBlock Chemi

- Protein-Free Blocking solution
- Enhanced Blocking efficiency than conventional product
- The largest volume among the series



EzBlock BSA

- The most widely usable Blocking solution including BSA
- Adjustable volume of Tween-20
- Consist of Tween-20 as a standard



EzBlock CAS

- The most effective protein based blocking solution including Casein
- Adjustable volume of Tween-20
- Consist of Tween-20 as a standard



Specification

Model	AE-1475	AE-1476	AE-1477
Product Name	<i>EzBlock Chemi</i>	<i>EzBlock BSA</i>	<i>EzBlock CAS</i>
Code	2332615	2332616	2332617
Package	1 Bottle	1 Set	1 Set
Volume	500mL(5-fold concentration) Mini size membrane 50pcs applicable	200mL(5-fold concentration) Tween-20 : 10mL Mini size membrane 20pcs applicable	200mL(5-fold concentration) Tween-20 : 10mL Mini size membrane 20pcs applicable
Components	TBS, Blocking reagent	TBS, BSA	TBS, Casein
Application	Western Blotting / ELISA / Anti-body dilution / Avidin-Biotin reaction / Analysis of Phosphoprotein / etc.	Western Blotting / ELISA / Anti-body dilution / Immunostaining / Avidin-Biotin reaction / Analysis of Phosphoprotein / etc.	Western Blotting / ELISA / Anti-body dilution / immunostaining / Southern Blot/ Northern Blot / etc.
Storage	Cold storage	Cold storage	Cold storage



ATTO Corporation

3-2-2 Moto-asakusa, Taito-ku, Tokyo 111-0041, JAPAN
Tel: +81-3-5827-4863 Fax: +81-3-5827-6647
e-mail: eig@atto.co.jp
URL: <http://www.atto.co.jp/eng/index.html>